

Jadelle®

Contraceptive Implants up to 5 years


First Implants WHO-prequalified*

*September 2009

Insertion




The implants will be inserted subdermally.

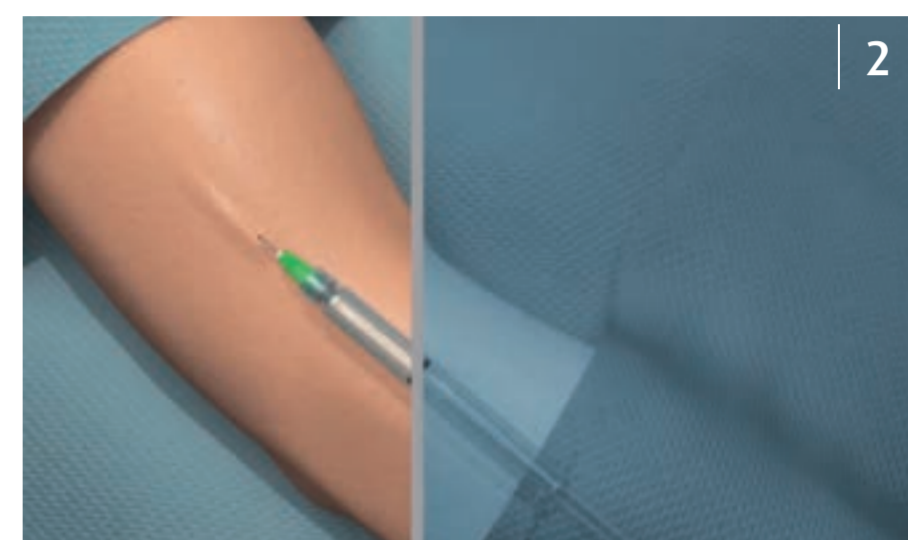


Anaesthetize the insertion area.




Use the new disposable trocar.

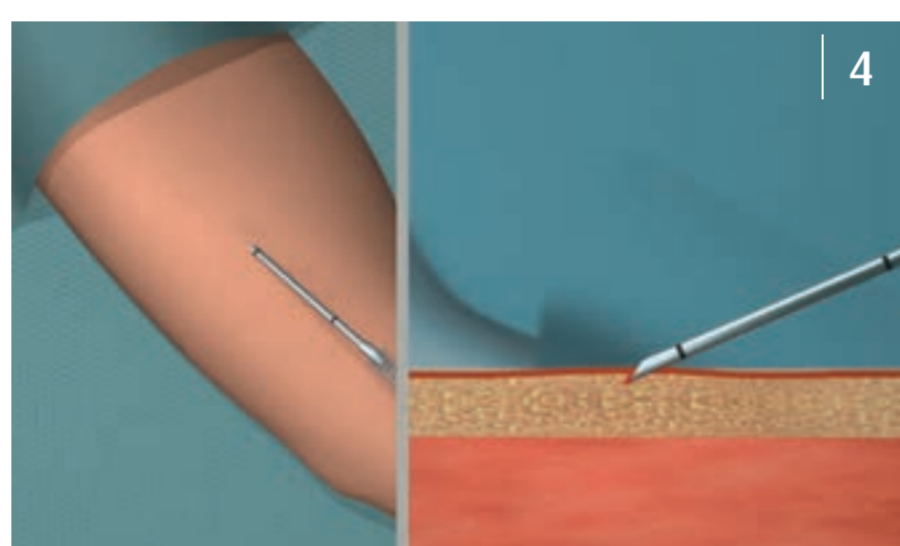


Introduce the trocar just beneath the skin.

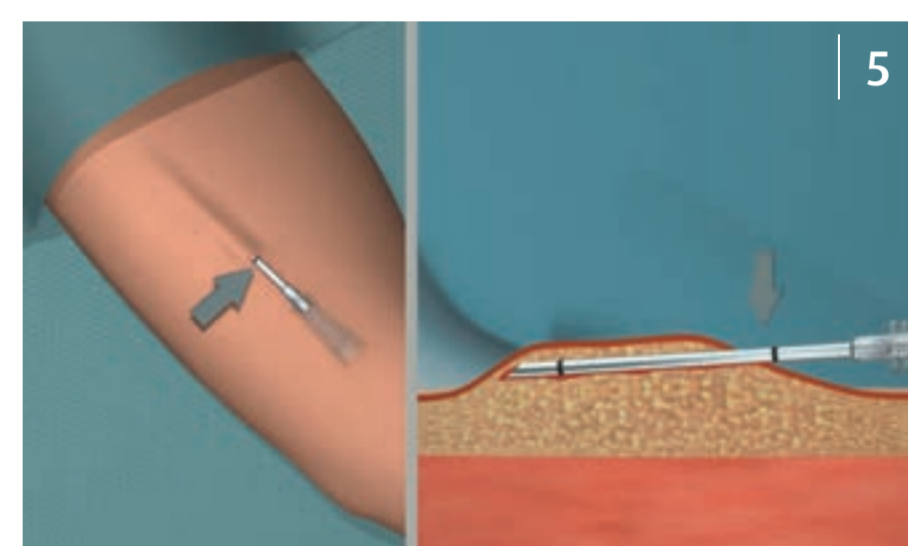


Advance to the ring mark, while "tenting".

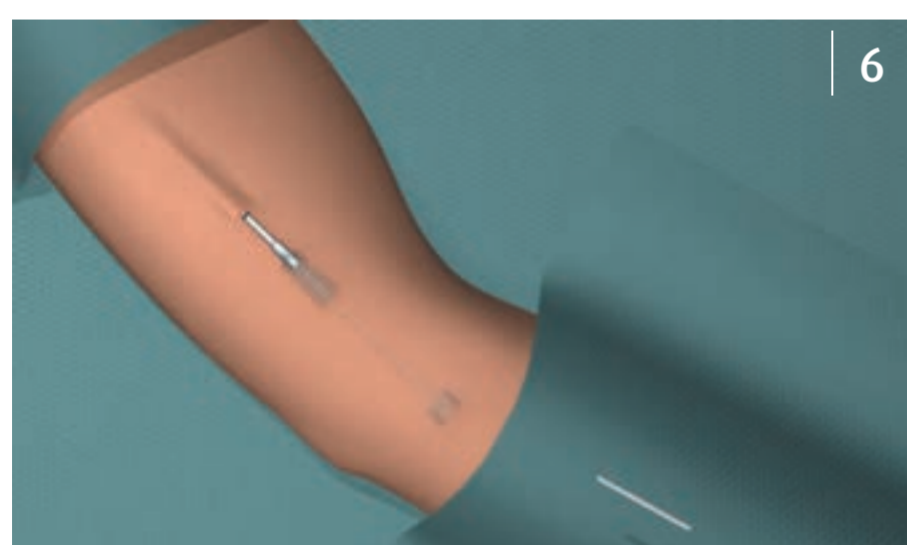


Remove the plunger and load the first implant.

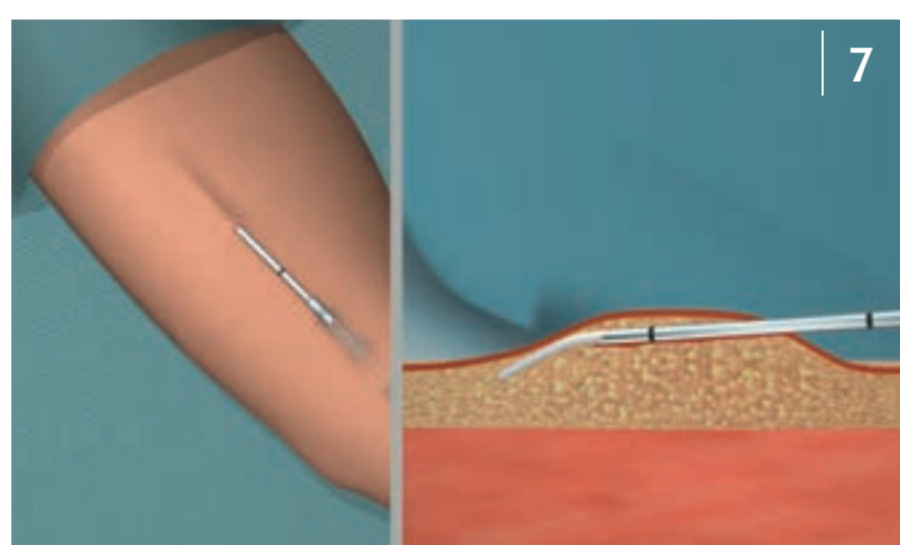


Hold the plunger steady and pull the trocar to the mark near the tip.

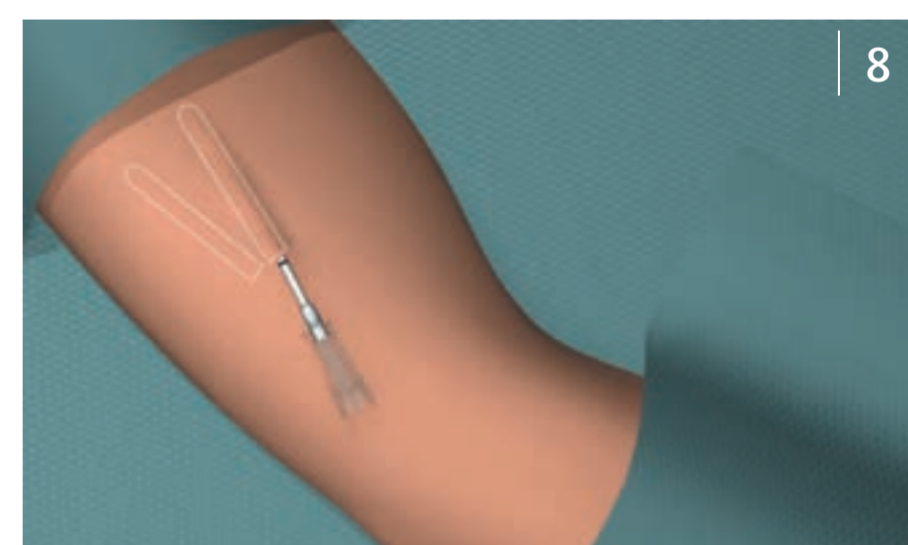


Advance again to the mark forming a narrow V.

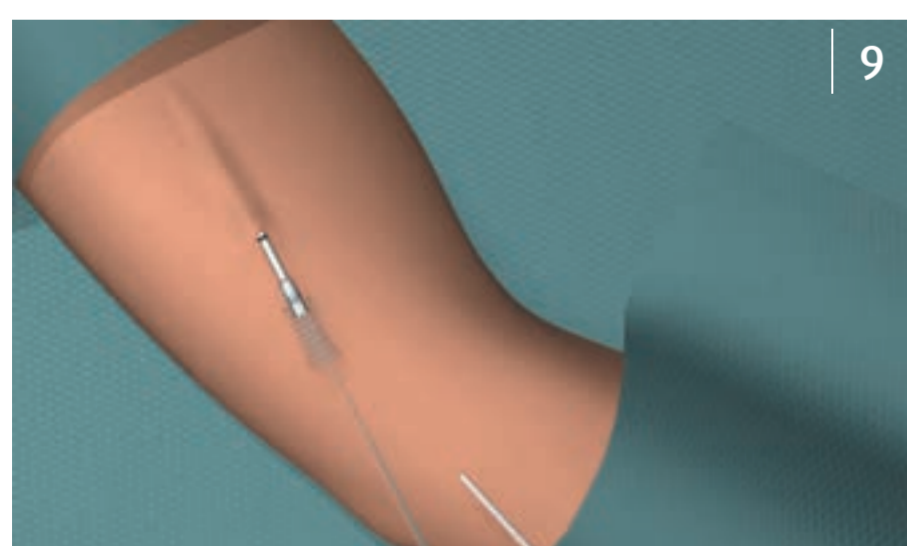

Remove the plunger and load the second implant.

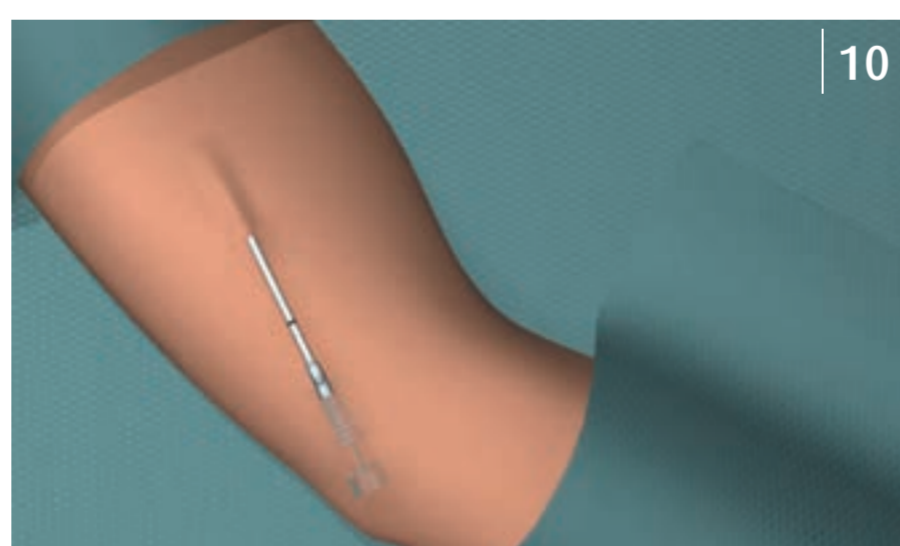

Hold the plunger steady while pulling the trocar back.

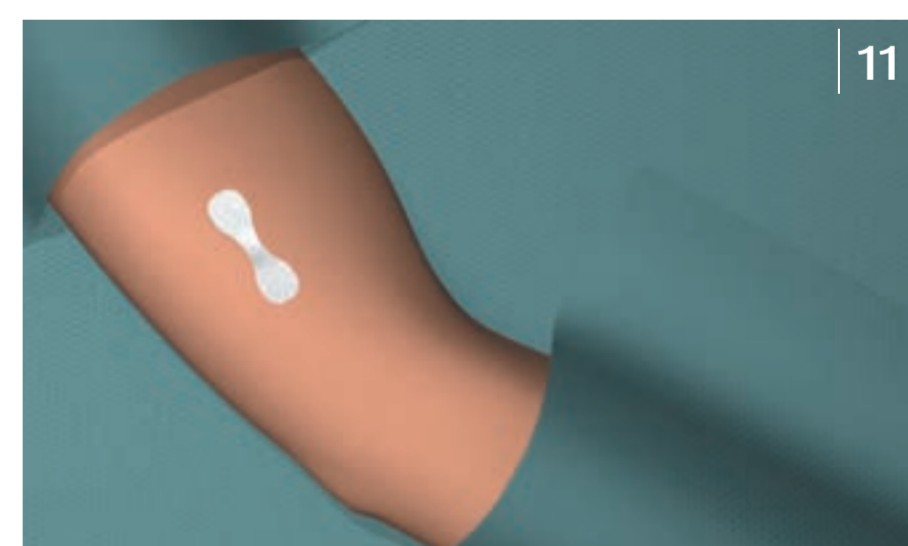

After the insertion close the incision with a sterile skin closure.